BENTLON

USER MANUAL



Bentlon Treatment bed SoSoft

60070 / 60071 (T) / 60078 / 60079 (T)

BENTLON

Treatment bed So Soft

CE

Bentlon / B Industries BV
De Noord 27
6001DA Weert / The Netherlands

Phone:	+31 (0)495 745 005
E-mail:	info@bentlon.com
Web:	www.bentlon.com

PREFACE

Dear customer,

We would like to take this opportunity to thank you for purchasing our quality product and for the trust you have placed in us.
Please read this instruction manual before using the So Soft for the first time. Reading the complete instruction manual is essential for the proper and safe use of the treatment bed SoSoft and will help you to operate the treatment bed SoSoft to your complete satisfaction.
Our treatment bed SoSoft is very robust and has a long lifespan when used properly.
Read the instructions carefully and operate the treatment bed SoSoft exclusively in accordance with the intended use as described in this instruction manual. This is the best to ensure safe usage of the treatment bed SoSoft and that improper use is avoided. Lack of knowledge on how to properly use the treatment bed SoSoft could result in problems for the customer and breakage of the treatment bed SoSoft.
Always keep the instruction manual within easy reach of the So Soft.

THIS INSTRUCTION MANUAL MUST ALWAYS BE ACCESSIBLE TO EVERY PERSON WHO OPERATES THE TREATMENT BED SOSOFT. MAKE SURE THAT THE INSTRUCTION MANUAL HAS BEEN READ AND UNDERSTOOD BEFORE OPERATING THE SPA TABLE.

⚠ WARNING! Non-compliance with the instructions in this manual endangers the operator's and customer's health and life and / or can damage the So Soft and will render the warranty invalid.

BENTLON

Index

1. Safety information 4
1.1 Abbreviations 4
1.2 Signal words 5
1.3 Safety instructions 5
2. Intended Use 8
3. Scope of Delivery 8
4. Preliminary remarks 8
5. Delivery and unpacking 9
5.1 Delivery 9
5.2 Unpacking 9
6. Assembly and transport 10
6.1 Assembly 10
6.2 Transport 11
6.2.1 Use of the carrying aids 12
7. Rating plate 14
8. Upholstery and design 14
8.1 Upholstery 14
8.2 Versions of the base 15
9. Operation 16
9.1 Characteristics 16
9.2 Model variants / Number of actuators 16
9.2.1 Treatment bed SoSoft V2 / 4 actuators 17
9.3 Before first use 17
9.4 View 17
9.5 Gas springs 18
9.6 Headrest 19
9.7 Armrests 20
10.7.1 Additional locking of the armrests 20
9.8 Legrest 21
9.9 Backrest 21
10. Electrical adjustments 21
10.1 Hand control 21

11.1.1 Adjustment of the working height ... 24
10.1.2 Adjustment of the backrest ... 24
10.1.3 Adjustment of the inclination of the seat section ... 25
10.1.4 Adjustment of the Trendelenburg / anti-Trendelenburg function ... 25
10.1.5 Memory functions of the hand control ... 25
10.1.6 Locking function of the hand control ... 26
10.2 Operation ... 28
11. Reset ... 29
12. Cleaning and Care ... 30
12.1 Cleaning of the upholstery ... 31
16.1.1 Cleaning and care of the polyurethane material (PU) ... 31
12.1.2 Cleaning and care of the artificial leather ... 31
12.2 Cleaning the housing and the fittings ... 31
12.3 Cover of the upholstery ... 31
13. Malfunctions (possible causes and troubleshooting) ... 32
13.1 Repair ... 33
14. Optionally available accessories ... 34
14.1 Protective covers made of terry ... 34
15. Upholstery care and cleaning agents ... 34
16. Conditions of transport, storage, and operation ... 34
16.1 Transport ... 34
16.2 Ambient conditions for storage and transport ... 34
16.3 Ambient conditions for operation ... 35
17. Technical Data ... 35
18. Disposal ... 36
19. Warranty ... 37
19.1 Exclusions from warranty ... 37
20. EC-Declaration of Conformity ... 38



Do not dispose of the treatment bed SoSoft with your domestic waste.
The treatment bed must be disposed of at a collecting centre for recyclable electric and electronic appliances. Do not remove the symbol / stickers from the treatment bed SoSoft.

1. Safety information

You must always observe basic safety precautions, including the following, when operating electrical products, especially when children are nearby:

BEFORE INITIAL OPERATION OF THE SO SOFT,
READ THE COMPLETE SAFETY INFORMATION
AND SAFETY INSTRUCTIONS ON THE FOLLOWING PAGES.

- **All information contained therein serve for the protection of the operator. If you ignore the safety information and safety instructions, you will endanger your health and life and as the health and life of your clients.**

MAKE SURE THAT THIS INSTRUCTION MANUAL
IS AVAILABLE TO ANYONE
OPERATING THE BENTLON TREATMENT BED SOSOFT
AND MAKE SURE THAT IT HAS BEEN READ
AND UNDERSTOOD
BEFORE OPERATING THE TREATMENT BED.

- This instruction manual is an integral part of the product and must be carefully kept in a safe place to consult it throughout the equipment service life.
  Copies of this instruction manual are available upon request.
- Bentlon / B Industries BV shall not be held responsible for any damages, direct or indirect, caused to people, animals or objects resulting from negligence or non-compliance with the instructions contained herein.
- Before disposing of a treatment bed, it is essential to **render it inoperable!** Always unplug the spa table from the mains and cut off the mains cable. Finally, you must dispose of it together with the mains plug.

★ **Keep this instruction manual at hand for when needed.**
Strictly observe the instructions to avoid damage to persons and property.

1.1 Abbreviations

SoSoft	➢ BENTLON SOSOFT BRONZE – SILVER – GOLD - PLATINUM
▪	➢ OPERATING ADVICE / TIP

1.2 Signal words

⚠ **DANGER!** indicates a hazardous situation which, if ignored, will result in death or serious injury.

⚠ **WARNING!** indicates a hazardous situation which, if ignored, could result in death or serious injury.

⚠ **CAUTION!** indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.

NOTICE! indicates a hazardous situation which, if not avoided, may result in damage to the treatment bed.

1.3 Safety instructions

⚠ DANGER

To reduce the risk of electrocution:

1. Non-compliance with the orders of this instruction manual endangers the operator's and customer's health and life and / or can damage the spa table.
2. All electrical work to connect the treatment bed SoSoft to the mains must be carried out by a qualified professional. Do not modify the energy supply. The connection must be carried out in accordance with the current local and legal regulations. RISK OF ELECTRIC SHOCK!
3. The mains voltage and the mains frequency of the mains supply **must** correspond to the values indicated on the rating plate.
4. Do **not** make any changes or manipulations to the Treatment bed SoSoft and/or its mains supply (mains cable, mains plug). In case of unauthorised changes or manipulations, **any** warranty of the B Industries B.V. is excluded.
5. **Do not connect** the treatment bed SoSoft to the mains if the Spa Table itself or the mains cable or the plug are visibly damaged. RISK OF ELECTRIC SHOCK!
6. The mains cable shall not be placed below the treatment bed SoSoft except it is expressly permitted by B Industries B.V.. RISK OF ELECTRIC SHOCK!
7. The mains cable must not be placed close to the actuators during connection or repositioning as it can be damaged while operating the lifting columns/actuators. RISK OF ELECTRIC SHOCK!
8. Never try to repair the treatment bed SoSoft yourself. If the treatment bed SoSoft does not operate properly, please contact our customer service or your dealer. Original spare parts should be used only.
9. If the mains cable and / or the plug are damaged, it must be replaced by B Industries B.V. or an authorized aftersales service or a qualified professional only. Never try to repair or replace the mains cable and / or the plug yourself. RISK OF ELECTRIC SHOCK!
10. **Always** unplug the treatment bed SoSoft if it is not used for a longer period.
11. Always disconnect the treatment bed SoSoft from the mains **before** cleaning and / or maintenance. RISK OF ELECTRIC SHOCK!
12. **Never** pull the mains cable to unplug the Treatment bed SoSoft. Always use the mains plug itself to unplug the Treatment bed SoSoft. RISK OF ELEC-TRIC SHOCK!
13. **Never** touch the mains plug, the mains switch, or other electrical components with wet or damp hands. RISK OF ELECTRIC SHOCK!
14. Do not place or store the treatment bed SoSoft close to a source of liquid such as a bath, a shower, or a washbasin. There is a danger that the treatment bed SoSoft may meet splashed water. RISK OF ELECTRIC SHOCK!


⚠ WARNING!

To reduce the risk of burns, electrocution, fire, or injury to persons:

1. Use a properly secured safety socket with the appropriate mains voltage and mains frequency. See technical data and rating plate.
2. **The electric components must not meet water.** RISK OF ELEC-TRIC SHOCK!
3. Check all cables regularly and contact our customer service if any cable is in improper condition. RISK OF ELECTRIC SHOCK!
4. The treatment bed SoSoft must not stand on the mains cable. Ensure that the mains cable is not damaged when moving the treatment bed SoSoft.
   Do not step on the mains cable. Do not stretch the mains cable.
5. Always ensure that the mains cable is not bent or squeezed during installation.
6. **Do not** use any extension cables, multiple sockets, or adapters. RISK OF FIRE!
7. **Do not** open the housings of the electric control panels! RISK OF ELECTRIC SHOCK!
8. Liquids **must not** get into the electric control panels! RISK OF ELECTRIC SHOCK!
9. Do not use water or any other liquids while carrying out any application. Do not spill water or any other liquids over, on or into the treatment bed SoSoft. Should you spill any water, **immediately** pull the mains plug and dry the Treatment bed SoSoft thoroughly. RISK OF ELEC-TRIC SHOCK!
10. **Never** immerse any parts of the treatment bed SoSoft in water.
    When cleaning the treatment bed SoSoft, make sure that no water penetrates the treatment bed SoSoft or its housing. RISK OF ELECTRIC SHOCK!
11. The cables of the models incl. the heating system and / or the Vibro Wellness 10 function[1] are between the upholstery and the under-frame. **The upholstery may only be removed by persons authorized by B Industries B.V..**
12. **For Treatment bed SoSoft with soft upholstery only**: for long periods of use of the heating system (optional equipment), always keep some free space between the back and seat upholstery. RISK OF BURNS DUE TO HEAT ACCUMULATION!



Correct

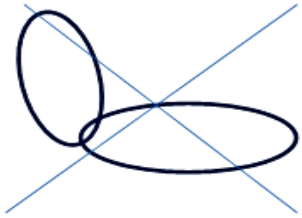

Wrong

13. **People with acute pain, venous disorders, low blood pressure, inflammatory processes and other medical conditions that make heat applications or this form of client positioning questionable must always consult a doctor before use.**
14. **Never use the Treatment bed SoSoft on people who are under the influence of alcohol or drugs.**
15. Operate or store the Treatment bed SoSoft in a dry environment **only**.
16. Close supervision is necessary if the treatment bed SoSoft is used by, on, or near children or invalids.
17. Use the treatment bed SoSoft only for the intended use as described in this instruction manual. Do not use attachments / accessories not recommended by B Industries B.V..
18. Never operate the treatment bed SoSoft if it has a damaged mains cable or plug, if it is not working properly, if it has been damaged or if it has become wet. Return the treatment bed SoSoft to an aftersales service authorized by B Industries B.V. for examination and repair.

19. In case of defects of the spa table, contact customer service **immediately**. Use of the spa table **must** be discontinued until the problem has been solved to avoid further damage to property and persons.
20. Regularly remove any dust build-up from the mains plug, the safety socket, and all electrical plug connections, e.g. with a commercially available vacuum cleaner. RISK OF FIRE!
21. Keep the mains cable away from heated surfaces.
22. Never drop or insert any object into any opening.
23. Do not use the Treatment bed SoSoft outdoors.
24. Do not allow any unauthorized person to operate the Treatment bed SoSoft!
25. Use the Treatment bed SoSoft only on a level surface.
26. **Do not place** any objects or liquids on the Treatment bed SoSoft.
    The Spa Table as well as the material can be damaged thereby.
    RISK OF ELECTRIC SHOCK! RISK OF DAMAGE!
27. Any repairs and maintenance work must be carried out by qualified professionals authorized by the B Industries B.V. **only**!
    RISK OF ELECTRIC SHOCK! RISK OF FIRE!
28. **Do not use** a steam cleaner for cleaning. RISK OF ELECTRIC SHOCK!
29. **Never** rest or stand on any part of the lying surface or of the supporting structure of the Treatment bed SoSoft. Do not exert strong physical force on these components; otherwise you can damage the supporting structure permanently. RISK OF TIPPING! RISK OF BREAKAGE! RISK OF INJURY!
30. Persons, animals or body parts of persons and animals **must not** be positioned between the base and the upper structure. RISK OF CRUSHING!
31. Ensure that no body parts or items are between the parts of the Treatment bed SoSoft while adjusting the lying position. RISK OF CRUSHING!
32. **Do not** use the supporting structure as a footrest. RISK OF CRUSHING!
33. Always ensure that the Treatment bed SoSoft stands firmly on its adjustable feet.
34. The Treatment bed SoSoft must **never** be used to carry or transport persons or heavy loads. RISK OF INJURY! RISK OF DAMAGE!
35. While unpacking, the packaging materials should be kept away from children, pets, or persons with reduced mental capabilities. RISK OF CHOKING! RISK OF INJURY!
36. Always keep children and pets away from the Treatment bed SoSoft.
37. Children must not play with the Treatment bed SoSoft or its single components.
38. The Treatment bed SoSoft is not intended to be operated by persons with reduced physical, sensory, and mental capabilities or lack of experience or knowledge, unless they are supervised or instructed concerning the use of the Treatment bed SoSoft and do understand the danger resulting thereof.

NOTICE!

1. The special pallet the Treatment bed SoSoft is delivered on has recesses. Upon delivery, the adjustable feet of the Treatment bed SoSoft stand in the respective recesses to prevent the Spa Table from inadvertent movement. Therefore, **it is essential to lift the Treatment bed SoSoft vertically off the special palle**t **using all four carrying aids**. Non-compliance with this instruction can seriously damage the adjustable feet of the Treatment bed SoSoft. RISK OF DAMAGE!
2. For safety reasons, the Treatment bed SoSoft only moves while keeping a button of the hand controlpressed. When released, the movement will stop immediately.
3. **Never** rest, stand, or kneel on any part of the lying surface of the Treatment bed SoSoft. Do not exert strong physical force on the lying surface. Otherwise you can damage the upholstery permanently. RISK OF DAMAGE!
4. Do not lean, stand, or sit on the armrests or headrest; otherwise you can damage the armrests or headrest permanently!
   RISK OF DAMAGE!
5. The maximum load capacity of the armrests is 25 kg.

6. To prevent unauthorized use of the Treatment bed SoSoft, **always** lock the drives manually using the hand control or the optional foot control when leaving the treatment room. No automatic locking is set as standard. You can decide for yourself whether to change this. Manual locking **does not** affect the heating cycle currently running and has **no effect** on the operation of the heating system.
   The heating cycle currently running is **not affected** by automatic locking.
7. The rating plate must not be removed or made illegible.
   Non-compliance will expire the warranty!

READ THIS INSTRUCTION MANUAL THOROUGHLY
AND STORE IT SAFELY.

2. Intended Use

Read the instructions very carefully and **only** operate the Spa Table in the intended use as described to protect yourself and your customers from danger.

- The Bentlon Treatment bed SoSoft is **no** medical device.
- The Treatment bed SoSoft is a versatile Spa Table. You can use it for e.g. cosmetic applications or full body massages in dry rooms, etc. **Any other use is not permitted and will affect the warranty claims.**
- The Treatment bed SoSoft is intended for use in closed rooms with defined ambient conditions. See chapter 20. CONDITIONS OF TRANSPORT, STORAGE, AND OPERATION.

★ **No liability for personal injury or material damage will be accepted if these instructions are not followed.**

3. Scope of Delivery

- 1 SPA TABLE Treatment bed SoSoft in desired design and equipment
- 1 HAND CONTROL
- 4 CARRYING AIDS
- 1 INSTRUCTION MANUAL

If you have purchased additional optional accessories for your Spa Table, these are also included in the scope of delivery.

4. Preliminary remarks

- The BentlonTreatment bed SoSoft is subject to the highest safety and quality standards.
- Before first use, make sure that all parts are present and that there is no visible damage.
- Make sure that the socket for the spa table is properly grounded and has the correct mains voltage and mains frequency, see rating plate.

.	**Before setting up, make sure that there is enough space for the Treatment bed SoSoft so that it does not bump into anything even when fully extended.**

5. Delivery and unpacking

5.1 Delivery

- The Treatment bed SoSoft is delivered and packed in a cardboard on a special pallet (s. fig. below). The Treatment bed SoSoft is in complete horizontal position and the upholstery is installed.



DELIVERY BY B INDUSTRIES B.V.:

If the Treatment bed SoSoft is delivered by the B Industries B.V. or one of our distributors, the delivery driver will help with the installation of the Treatment bed SoSoft and will also dispose of the packaging material and the special pallet.


DELIVERY BY AN EXTERNAL FORWARDER:

If the Treatment bed SoSoft is delivered by an external forwarder, the Treatment bed SoSoft is usually **only** delivered (no help with the installation of the Treatment bed SoSoft).
In this case, you must dispose of the packaging material and the special pallet on your own responsibility.

5.2 Unpacking

⚠ **WARNING!** While unpacking, the packaging should be kept away from children, pets, or persons with reduced mental capabilities. RISK OF CHOKING! RISK OF INJURY!

- The Treatment bed SoSoft is delivered and packed in a cardboard on a special pallet (s. fig. below).

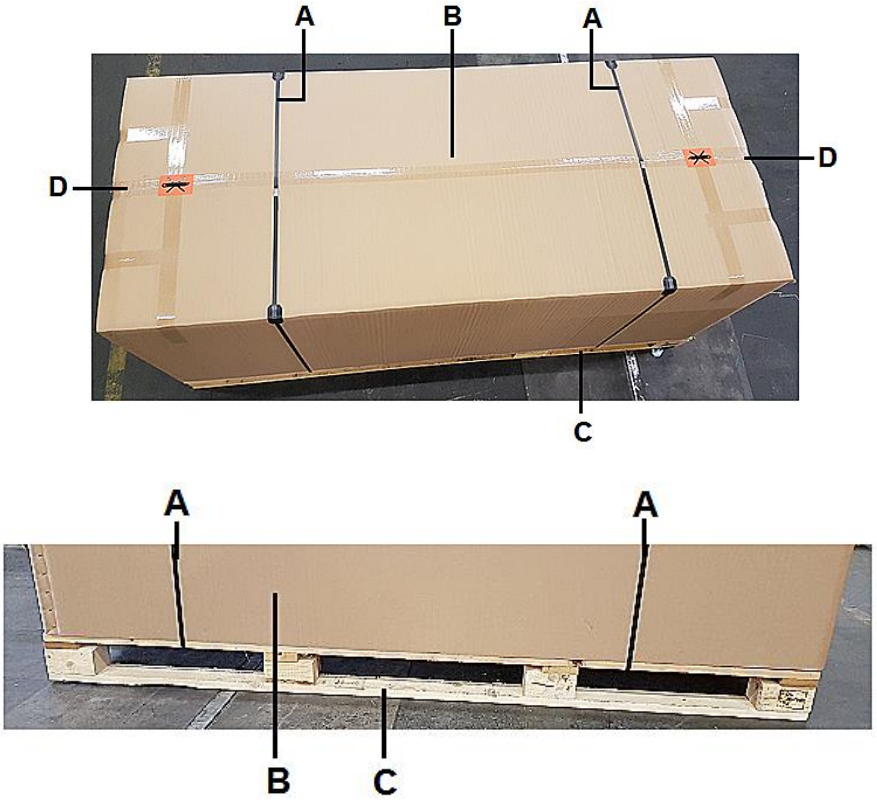

A	PACKAGING TAPES
B	UPPER CARDBOARD
C	SPECIAL PALLET
D	SHORT SIDES OF THE CARDBOARD

1. Cut the packaging tapes (**A**) with scissors and remove them.
2. **DO NOT open** the upper cardboard (**B**) **but** lift it carefully and vertically as a whole from the Treatment bed SoSoft and from the special pallet (**C**). For this purpose, you will need two persons; one person each on the short sides (**D**) of the cardboard.
3. Dispose of the cardboard properly.
4. Now you can move the Treatment bed SoSoft from the special pallet; **strictly observe** the instructions given in chapter 6.2 TRANSPORT / CARRYING AIDS and chapter 6.2.1 USE OF THE CARRYING AIDS when doing so.

6. Assembly and transport

6.1 Assembly

The Treatment bed SoSoft is usually delivered fully assembled, so no assembly work is necessary. In **exceptional cases**, necessary assembly work will be carried out by the delivery driver on site.

The Treatment bed SoSoft is easy to install and can be adapted to all human body measurements. You can adjust the desired positions electrically in a very short period and then focus entirely on your customer.

- Before setting up, make sure that there is enough space for the Treatment bed SoSoft so that it does not bump into anything even when fully extended.

⚠ **WARNING!** The cables of the models incl. the heating system is between the upholstery and the underframe.**The upholstery may only be removed by persons authorized by B Industries B.V. or distributors of B Industries B.V.!**

6.2 Transport

- Transport of the Treatment bed SoSoft can be severely affected to some extent due to unfavourable environmental conditions. Such unfavourable environmental conditions are e.g. narrow doors, narrow corners behind doors, narrow corridors, narrow staircases, spiral staircases etc. The weight of the Treatment bed SoSoft is approx. 120 kg to 150 kg (depending on the design of the base); the length is approx. 200 cm and the width including the upholstery is 75,4 cm / 80 cm/ 85 cm or 91,5 cm^2; thus, the manoeuvrability of the Treatment bed SoSoft is limited

CARRYING AIDS:

- You need the carrying aids (s. fig. below) **to lift the Treatment bed SoSoft from the special pallet**. The proper use of the carrying aids is described below in chapter 6.2.1.
  You can optionally purchase the carrying aids as a set.

- **You need 4 persons** to move or transport the Treatment bed SoSoft using the carrying aids.

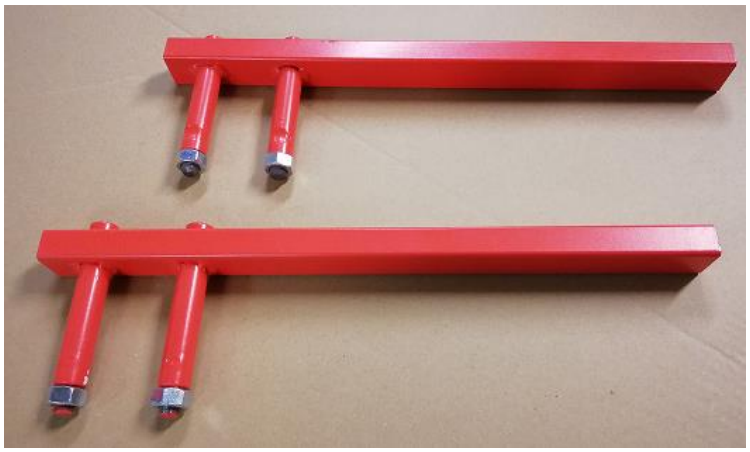


Fig. similar: fig shows 2 of a total of 4 carrying aids.

6.2.1 Use of the carrying aids

⚠ DANGER! After transport, while putting down the Treatment bed SoSoft, the mains cable must not be under the Spa Table; otherwise the adjustable feet may squeeze the mains cable and cause permanent damage while putting down the Treatment bed SoSoft. RISK OF ELECTRIC SHOCK!

NOTICE! The special pallet the Treatment bed SoSoft is delivered on has recesses (s. figure). Upon delivery, the adjustable feet of the Treatment bed SoSoft stand in the respective recesses to prevent the Spa Table from inadvertent movement. **Therefore, it is essential to lift the Treatment bed SoSoft vertically off the special pallet** using all four carrying aids. Non-compliance can seriously damage the adjustable feet of the Treatment bed SoSoft. RISK OF DAMAGE!




Figures similar: all following figures serve as an example only.


1. To move the Treatment bed SoSoft from the special pallet, you must raise the Treatment bed SoSoft first. Plug in the Treatment bed SoSoft to do this. Then press the corresponding button of the hand control which is on one side on a suspension device / (**button T 1**, s. fig. in chapter 11.1 HAND CONTROL.

2. You must bring the **backrest** into an upright position before you use the carrying aids. To do this, press the corresponding button of the hand control (**button T 2**, s. fig. in chapter 11.1 s. fig. in chapter 11.1 HAND CONTROL.

3. **Then unplug the Treatment bed SoSoft.**

4. Insert the 4 carrying aids (**A**) into the intended openings (**B** / 2 front left, 2 front right, 2 back left, 2 back right) underneath the lying surface (**C**); s. fig. below.

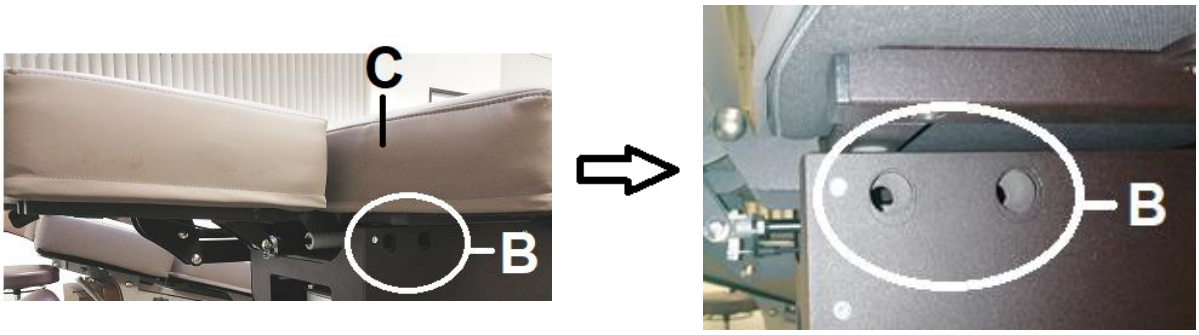

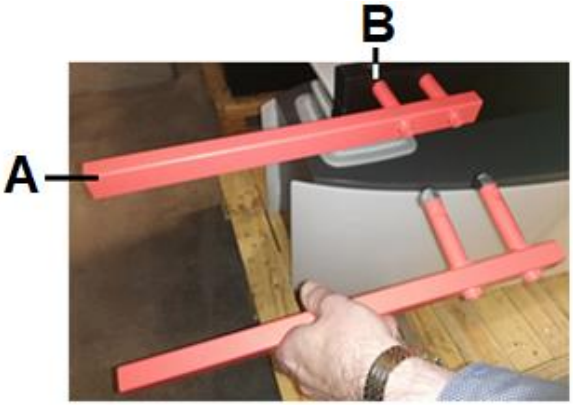



A	CARRYING AIDS
B	OPENINGS FOR THE CARRYING AIDS
C	LYING SURFACE

5. You need **4 persons** to move the Treatment bed SoSoft to its final site.

6. It is essential to **lift** the Treatment bed SoSoft **vertically** off the special pallet using **all four** carrying aids. Otherwise you damage the adjustable feet.

7. Carry or transport the Treatment bed SoSoft **using all four carrying aids only**; otherwise you can damage the Treatment bed SoSoft irreparably.

8. Store the carrying aids in a safe place. If the Treatment bed SoSoft must be carried or transported at a later point of time (e.g. when moving), it is only possible when using the four carrying aids.

9. If the floor is not level, you must ensure proper stability of the Treatment bed SoSoft by adjusting the four adjustable feet (**D** /s. fig. below). After you have adjusted the feet, they must be secured by tightening the respective lock nut (**E**); see fig. below.



D	ADJUSTABLE FOOT
E	LOCK NUT

NOTICE! Once you have found a suitable place with enough space for the spa table, remove the carrying aids and store them well for possible later transport!

NOTICE! If the floor is not level, you **must** ensure proper stability of the Treatment bed SoSoft by adjusting the four adjustable feet.


7. Rating plate

- ❖ The rating plate is located at the head end of the Treatment bed SoSoft; on the wooden plate below the backrest (see figure below).
- ➢ To gain access to the rating plate, you must bring the backrest of the Treatment bed SoSoft into an upright position.

Do not make the rating plate illegible! Do not remove the rating plate! Noncompliance will expire the guarantee!

·	**If you wish to contact the Bentlon customer service regarding your Spa Table, we recommend taking a photo of the rating plate previously (e.g. using a mobile phone or a digital camera). This makes it easier for you to identify the information of the rating plate (e.g. the REF and CE number of the Spa Table) to provide the required information to the customer service. On the cover page of this instruction manual is an area where you can note the serial number of your Spa Table. This way you always have the serial number at hand.**

1. You find the REF en CE information on the silver Bentlon label. The label sticks to one of the actuatorkoloms.
2. Below you find an example of the silver Bentlon Label.

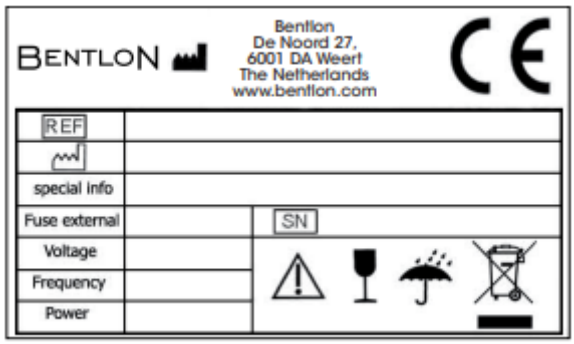

8. Upholstery and design

8.1 Upholstery

- ❖ The Treatment bed SoSoft is available in **one** upholstery variants
- ❖ The upholstery of the Treatment bed SoSoft is made of **polyurethane** (**PU**).

STANDARD UPHOLSTERY MADE OF POLYURETHANE (PU):

- The Treatment bed SoSoft can be equipped with a soft, comfortable, and high-quality multi-piece upholstery (s. fig. below).
- The inner material of the upholstery is made of high-quality foam material.
- The upholstery is made of **polyurethane** (**PU**) and available in various colours.



Fig. similar: colour and decor can vary depending on model.

8.2 Versions of the base

- The Treatment bed SoSoft is available with two different version of the base.

TREATMENT BED SOSOFT BRONZE AND GOLD:

- The base of the Treatment bed SoSoft Bronze and Gold consists of a solid plate with elegant wooden decor (s. fig. below).
- The wood base is available in five different colours.



9. Operation

⚠ **WARNING!** Do not sit on the headrest and armrests.

⚠ **WARNING!** Do not, under any circumstances, exert physical force on the above-mentioned components, regardless of the lying surface setting. Non-compliance with this instruction may permanently damage the support structure of the Treatment bed SoSoft. RISK OF TIPPING! RISK OF BREAKAGE! RISK OF INJURY!

9.1 Characteristics

- The Treatment bed SoSoft is a versatile Spa Table. You can use it for e.g. cosmetic applications, podiatry applications or full body massages in dry rooms.
- The Treatment bed SoSoft is available with two different versions of the base (TREATMENT BED SOSOFT BRONZE AND GOLD / TREATMENT BED SOSOFT SILVER AND PLATINUM).
- The Treatment bed SoSoft is available in two upholstery variants (s. chapter 9.1 UPHOLSTERY).
- The entry height of the Treatment bed SoSoft (**soft upholstery**) is approx. **+ 5 cm** compared to the other upholstery variants due to the thicker upholstery.
- The Bentlon Treatment bed SoSoft is a Spa Table with four actuators, which allow you to adjust the spa table with a hand control.

- Frequently used treatment positions can be adjusted, stored, and recalled by two individual memory locations (s. chapter 10.2 MODEL VARIANTS / NUMBER OF ACTUATORS) using the hand control (s. chapter 11.1.5 MEMORY FUNCTIONS OF THE HAND CONTROL).
- Working height, backrest, inclination of the seat section and the Trendelenburg / Anti-Trendelenburg function are electrically adjustable, see chapter 11. ELECTRICAL ADJUSTMENTS. The inclination of the seat section is adjustable by the electrically adjustable knee bend.
- The armrests can be adjusted pneumatically; s. chapter 10.7 ARMRESTS.
- The headrest can be adjusted pneumatically; s. chapter 10.6 HEADREST.
- The maximum lifting capacity of the Treatment bed SoSoft is approx. 300 kg.
- The maximum load capacity of the lying surface of the Treatment bed SoSoft is approx. 350 kg.
- The maximum load capacity of the armrests is 25 kg.

9.2 Model variants / Number of actuators

- The Spa Tables of the Bentlon Treatment bed SoSoft series are available in one size widths: SMALL 78 CM.

- The Spa Tables of the Bentlon Treatment bed SoSoft series can be equipped with four actuators, which allow you to adjust the spa table with a hand control.

9.2.1 Treatment bed SoSoft V2 / 4 actuators

Treatment bed SoSoft V2/ 4 ACTUATORS:	
2 lifting columns	Electrical adjustment of the WORKING HEIGHT. Electrical adjustment of the TRENDELENBURG / ANTI-TRENDELENBURG FUNCTION
Actuator 1	Electrical adjustment of the BACKREST.
Actuator 2	Electrical adjustment of the INCLINATION OF THE SEAT SECTION: adjustable by the electrically adjustable KNEE BEND.

- The LEGREST is **not** adjustable.

9.3 Before first use

- Before first use, please check that all components are present (s. chapter 3. SCOPE OF DELIVERY) and that they are not visibly damaged.
- Make sure that the socket for the spa table is properly grounded and has the correct mains voltage and mains frequency, see rating plate.

9.4 View

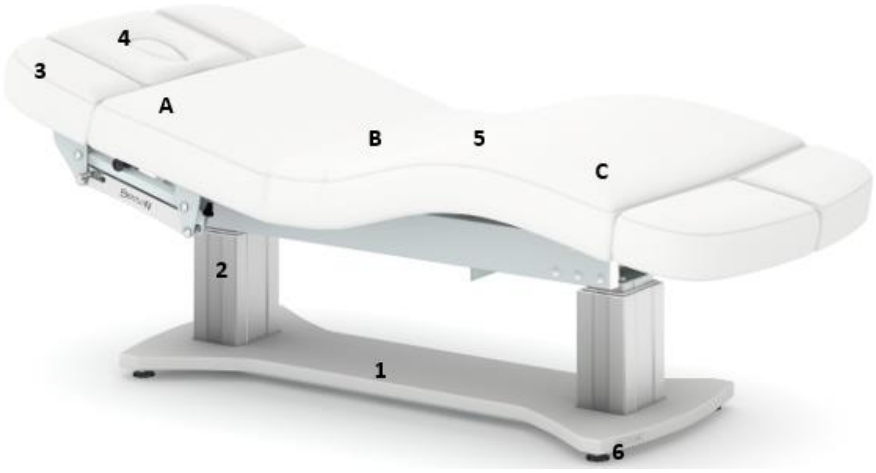


Fig. similar: colour, design and equipment may vary depending on model.

1	**BASE**

2	**HEIGHT ADJUSTMENT / INCLINATION ADJUSTMENT:** ➢ s. chapter 11.1.1 ADJUSTMENT OF THE WORKING HEIGHT ➢ s. chapter 11.1.4 ADJUSTMENT OF THE TRENDELENBURG / ANTI-TRENDELENBURG FUNCTION
3	**ARMRESTS:** ➢ Pneumatically adjustable and swivelling. ➢ s. chapter 10.5 GAS SPRINGS and chapter 10.7 ARMRESTS
4	**HEADREST:** ➢ Pneumatically adjustable. ➢ Incl. nose slot. ➢ s. chapter 10.5 GAS SPRINGS and chapter 10.6 HEADREST
5	**LYING SURFACE, DIVIDED INTO:** **a. Electrically adjustable backrest** ➢ s. chapter 10.9 BACKREST and chapter 11.1.2 OF THE BACKREST **b. Knee bend (inclination of the seat section)** ➢ The INCLINATION of the seat section is adjustable by the electrically adjustable KNEE BEND. ➢ s. chapter 11.1.3 ADJUSTMENT OF THE INCLINATION OF THE SEAT SECTION **c. Legrest** ➢ The legrest **is not** adjustable.
6	**ADJUSTABLE FEET:** ➢ s. chapter 6. ASSEMBLY AND TRANSPORT

9.5 Gas springs

❖ The armrests and the headrest can be adjusted pneumatically with a gas spring using a lever.

❖ If you carry out any setting using the lever of the gas spring, you **only** need to operate the lever **when lowering** the armrests or headrest.

- Carry out a pulling movement upwards (**without operation of the respective leve**r) **to lift** an armrest and / or the headrest again.

9.6 Headrest

⚠ **WARNING!** Do not carry out any physical force to the headrest, regardless of the setting of the lying surface. Non-compliance can cause permanent damage to the construction of the Treatment bed SoSoft. RISK OF TIPPING! RISK OF DISINTEGRATION! RISK OF INJURY!

The headrest is freely and manually adjustable (s. below: **fig. 1, 2 and 3**).

You can incline the headrest (positive by 35° and negative by 48°) using the lever, which is on the underside of the headrest (**fig. 1**).
Lock the set position by releasing the lever.

- see also chapter 10.5 GAS SPRINGS



FIG. 1



FIG. 2



FIG. 3

Fig similar: figures serve as an example only.

- The upholstery of the headrest is equipped with a nose slot usable when carrying out back massages. For this purpose, the headrest is equipped with
a separate upholstery part, which can be removed (s. fig. 4).

9.7 Armrests

⚠ **WARNING!** Do not carry out any physical force to the armrests, regardless of the setting of the lying surface. Non-compliance can cause permanent damage to the construction of the Treatment bed SoSoft. RISK OF TIPPING! RISK OF DISINTEGRATION! RISK OF INJURY!

- You can lower and lift each armrest in the same way as the headrest (manually, gas spring / see also chapter 10.5 GAS SPRINGS).
- Additionally, the armrests of the Treatment bed SoSoft can also be swivelled and locked in at 45° and at 90°.
- Use the lever of the respective side and press lightly to adjust the armrests; so the Treatment bed SoSoft is the optimal solution for back massages (s. **fig. 10**).

As you can swivel the armrests, you are closer to the head of the person to be treated when carrying out cosmetic or facial massages.



FIG. 10

Fig similar: figures serve as an example only.

NOTICE! Swivel the armrests without load only! Do not lean, stand, knee, or sit on the armrests. Otherwise you can damage the armrests permanently. RISK OF DAMAGE!

NOTICE! The maximum load capacity of the armrests is 25 kg. RISK OF DAMAGE!

10.7.1 Additional locking of the armrests

- The armrests of the Treatment bed SoSoft can be locked **additionally** using the laterally mounted locking button (s. fig. below).
- The additional lock enables higher stability and safety of the armrests.

LOCKING OF THE ARMRESTS:

1. **Push** the inner **red** button **inward** to lock the armrests.
2. Now the armrests are locked.

UNLOCKING OF THE ARMRESTS:

1. **Pull** the outer **black** button **forward** (towards you) to unlock the armrests.
2. The additional lock of the armrests is released.



NOTICE! The additional lock must be released to adjust the armrests using the respective levers (s. chapter 10.7 ARMRESTS). Otherwise the construction and the mechanics of the armrests can be damaged permanently. RISK OF DAMAGE!

9.8 Legrest

The two single legrests are adjustable.

9.9 Backrest

- The backrest can be adjusted electrically, s. chapter 11.1.2 ADJUSTMENT OF THE BACKREST.

10. Electrical adjustments

- The electrical adjustments of the Treatment bed SoSoft are carried out with a hand control.
- The pneumatic adjustments of the Treatment bed SoSoft (armrests / headrest) are carried out with a gas spring using a lever.

·	**For safety reasons, the Treatment bed SoSoft only moves while keeping a button of the hand control pressed. When released, the move-ment will stop immediately.**

10.1 Hand control

The hand control is a wired control.
It enables operation of the electrically controllable functions.


- The two lifting columns and the actuators can be controlled by pressing and holding down the corresponding button.
- Furthermore, two memory positions can be adjusted, stored and recalled. How to save certain lying positions is described in chapter 11.1.5 MEMORY FUNCTIONS OF THE HAND CONTROL.
- Recalling the four preset positions is also possible.
- The hand control has 12 buttons.
- Adjustment of the working height, the backrest, the inclination of the seat section, and the Trendelenburg / anti-Trendelenburg function (buttons T 1 - T 4) is done by a double-tap

system. This is for safety reasons.
A movement that is associated with danger is made consciously by this.

For the direction of movement, please note the arrow symbols and the number of dots in the upper area of the hand control.

By simply pressing a button and keeping it pressed, a movement is made upwards.
By pressing a button twice (double-tap system) and keeping it pressed,
a movement is made in the opposite direction.

The two individual memory positions and de-activation of the locking function are also done via the double-tap system.

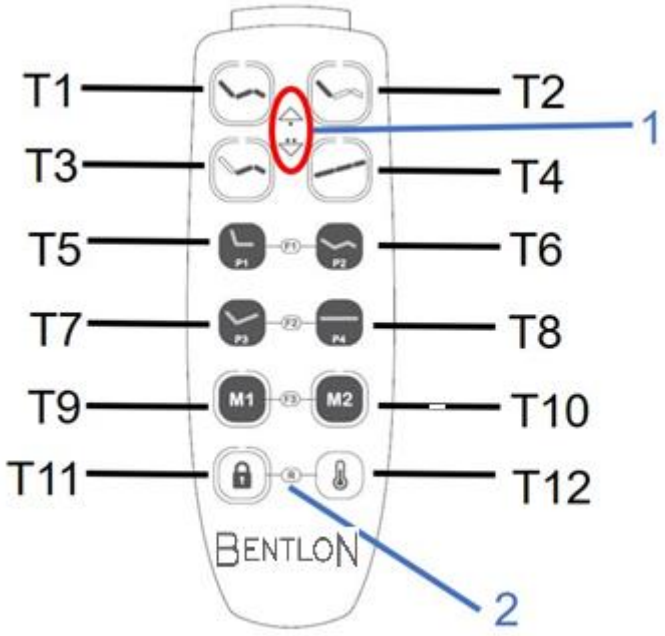

T 1	Button for adjustment of the **WORKING HEIGHT**
T 2	Button for adjustment of the **BACKREST**
T 3	Button for adjustment of the **INCLINATION** (seat area)
T 4	Button for adjustment of the **TRENDELENBURG / ANTI-TRENDELENBURG FUNCTION**
T 5	**PROGRAM BUTTON P1** ➢ PROGRAM BUTTON P1 moves the Treatment bed SoSoft to the fixed entry position. ➢ The fixed entry position cannot be changed.

T 6	**PROGRAM BUTTON P2** ➢ PROGRAM BUTTON P2 moves the Treatment bed SoSoft to the fixed S-curved position. ➢ The S-curved position cannot be changed.
T 7	**PROGRAM BUTTON P3** ➢ PROGRAM BUTTON P3 moves the Treatment bed SoSoft to the fixed pedicure position. ➢ The pedicure position cannot be changed.
T 8	**PROGRAM BUTTON P4** ➢ PROGRAM BUTTON P4 moves the Treatment bed SoSoft to the fixed massage position. ➢ The massage position cannot be changed.
T 9	**MEMORY BUTTON M1** ➢ MEMORY BUTTON M1 is used to save and recall the first individually adjusted and saved position.
T 10	**MEMORY BUTTON M2** ➢ MEMORY BUTTON M2 is used to save and recall the second individually adjusted and saved position.
T 11	**LOCK BUTTON** ➢ The LOCK BUTTON locks and unlocks the hand control / the actuators.
T 12	**BUTTON FOR THE HEATING SYSTEM** (optional) ➢ With this button, you can switch on the heating system, adjust it in 3 levels and switch it off again.

1	**SYMBOLS FOR THE DOUBLE-TAP SYSTEM** ❖ The working height, backrest, inclination of the seat section and Trendelenburg / anti-Trendelenburg function (buttons T 1 - T 4) are adjusted using the double-tap system. ❖ This serves for safety. A movement associated with danger is made consciously by this.
2	**SYMBOL FOR RESET** ❖ Press buttons T 11 + T 12 on the hand control simultaneously and keep them pressed. ❖ See chapter 15. RESET. ❖ **A reset procedure is only possible in the event of a system error.**

11.1.1 Adjustment of the working height

- For adjustment of the working height press the corresponding button (**T 1**); s. fig. in chapter 11.1 HAND CONTROL.
- For the direction of movement, please note the arrow symbols and the number of dots above them on the hand control.
- Press the button once and hold it down to raise the Treatment bed SoSoft. Press the button twice and hold it down to lower the Treatment bed SoSoft.
- The horizontal WORKING HEIGHT of the Treatment bed SoSoft is adjustable from approx. 66 cm up to approx. 96 cm, due to the thicker upholstery (if the adjustable feet are retracted).
- The entry height is the lowest position of the seat / seat section. At this height, your customer has easy access to the Treatment bed SoSoft from the side.
- Treatment bed SoSoft = approx. 66 cm / seat (seat section) in horizontal position

10.1.2 Adjustment of the backrest

- The inclination of the backrest can also be adjusted using the hand control.
- For adjustment of the backrest press the corresponding button (**T 2**); s. fig. in chapter 11.1 HAND CONTROL).
- For the direction of movement, please note the arrow symbols and the number of dots above them on the hand control.

- Press the button once and hold it down to raise the backrest of the Treatment bed SoSoft. Press the button twice and hold it down to lower the backrest of the Treatment bed SoSoft.

⚠ **WARNING!** Ensure that no body parts or items are between the parts of the Treatment bed SoSoft while adjusting the treatment position. RISK OF CRUSHING!

10.1.3 Adjustment of the inclination of the seat section

- Via the electrically adjustable knee bend the inclination of the seat section (thigh section to knee bend section) can also be adjusted using the hand control.
- For adjustment of the inclination of the seat section press the corresponding button (**T 3)**; s. fig. in chapter 11.1 HAND CONTROL.
- For the direction of movement, please note the arrow symbols and the number of dots above them on the hand control.
- Press the button once to move up or twice to move down and hold it down to adjust the inclination of the seat section.

10.1.4 Adjustment of the Trendelenburg / anti-Trendelenburg function

- For adjustment of the inclination of the lying surface press the corresponding button (**T 4**, s. fig. in chapter 11.1 HAND CONTROL).
- For the direction of movement, please note the arrow symbols and the number of dots above them on the hand control.

- The Treatment bed SoSoft can be inclined to the desired direction of movement up to an inclination angle of approx. 14 °.
- If an insufficient height is reached during the adjustment of the inclination, the Treatment bed SoSoft moves up automatically to avoid collision with the base.
- During the crossover from Trendelenburg to Anti-Trendelenburg, the Treatment bed SoSoft stops for approx. one second when reaching the horizontal position. If you continue to press the corresponding button, the Treatment bed SoSoft will continue to move.


⚠ **WARNUNG!** Ensure that there are no persons, body parts of persons, animals or objects between the lying surface and the base. RISK OF COLLISION! RISK OF CRUSHING! RISK OF INJURY! RISK OF DAMAGE!

10.1.5 Memory functions of the hand control

- In addition to four pre-set treatment positions, the hand control of the Treatment bed SoSoft has two individual memory positions (memory positions: M 1 / M 2).
- The memory function saves the treatment positions you wish to use and simplifies your work processes. Use the buttons M 1 and M 2 (s. fig. in chapter 11.1 HAND CONTROL).
- The memory positions are saved in the same way as the electrical settings of the Treatment bed SoSoft using the double-tap system.

SAVING THE MEMORY POSITIONS:	
1. Adjust the desired position as usual with the buttons T 1 to T 4.	

2. Press button M 1 or M 2 twice /(double-tap system) and keep it pressed for about 3 seconds.

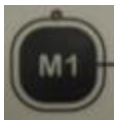

3. An acoustic signal sounds when the setting has been successfully saved.



4. The respective memory position is now stored.

RECALLING THE MEMORY POSITIONS:

Press and hold the corresponding memory button M 1 or M 2 pressed **until** the stored treatment position is reached.





10.1.6 Locking function of the hand control

- The Treatment bed SoSoft is equipped with a LOCKING FUNCTION to protect the spa table from unauthorized use.
- By default, no automatic locking is set. You can decide for yourself whether to change this.
- The currently running heating cycle is **not affected** by an automatic locking.
  A manual lock also **does not affect** the currently running heating cycle and has **no effect on the operation of the heating system.**

VARIOUS AUTOMATIC LOCKING TIMES:

Various automatic locking times can be set.
To get feedback on the currently set automatic locking time, press buttons T 9 and T 10 once simultaneously.

Feedback on the currently set automatic locking time is given as follows via flashing of all LEDs of the hand control ana several beeps, depending on the automatic locking time:

- automatic locking after **30** minutes:
  1 beep and **1 time** flashing of all LEDs of the hand control
- automatic locking after **60** minutes:
  2 beeps and **2 times** flashing of all LEDs of the hand control
- automatic locking after **90** minutes:
  3 beeps and **3 times** flashing of all LEDs of the hand control
- automatic locking after **120** minutes:
  4 beeps and **4 times** flashing of all LEDs of the hand control

- **no** automatic locking:
  5 beeps and **5 times** flashing of all LEDs of the hand control


To change from one automatic locking time to the next, please press both buttons T 9 and T 10 simultaneously and keep them pressed for about 6 seconds.
All LEDs flash and a corresponding beep is sounded to provide feedback on the currently set automatic locking time, followed by a change to the next automatic locking time.

Example:

30 minutes is set as the automatic locking time.
You want to change to the automatic locking time 60 minutes.
Press the corresponding buttons simultaneously and keep them pressed for approx. 6 seconds until feedback about the now set automatic locking time 60 minutes
is given through flashing of all LEDs of the hand control and beeps.

The automatic locking time is always changed in this sequence:

30 min. - 60 min. - 90 min. - 120 min. - **no automatic locking.**

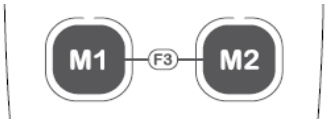

Buttons T 9 + T 10 (F3)

LOCKING THE DRIVES MANUALLY:

1. To lock the drives manually, press the LOCK button (**T 11**) and **keep it pressed** for about 2 seconds.
2. Successful locking is confirmed by an acoustic signal and the LED above the button on the hand control also lights up.
3. The drives are now locked.



UNLOCKING THE DRIVES MANUALLY:

1. To unlock the drives, press the LOCK button (**T 11**) on the hand control **twice in quick succession** (double-tap) and keep it pressed for about 2 seconds.

2. The successful unlocking is confirmed by an acoustic signal and the LED above the button on the hand control is switched off.

3. The drives are now **un**locked.



10.2 Operation

❖ The heating system is operated using button **T 12** of the hand control.



POWER ON:

❖ With this button to switch **on** the heating system.

The heating system comes with three heat / power levels.
To quickly heat up the lying surface before starting treatment, we recommend setting the heating system to level 2; only in exceptional cases to level 3 (highest level).

Please refer to the following charts to find information on how to operate the heating system.

After a total time of **90 minutes** the heating system switches off **automatically**; regardless of the set heating level.

❖ This process is pre-programmed.



POWER OFF:

❖ With this button you also switch **off** the heating system.

There are three heat / power levels to choose from:

		Heating level	Visualization via LED (diode)	Heating output
button	Push button **once**	Power on / level **1**	**single** impulse for 0.3 seconds, interval 1 second	30 %
button	Push this button **again**	level **2**	**double** impulse, lights up for 0.3 seconds, interval 1 second	60 %
button	Push this button **again**	level **3**	**triple** impulse, lights up for 0.3 seconds, interval 1 second	90 %
button	Push this button **again**	Power off		

NOTICE! If the heating system is not present on your spa table, this will be indicated by the LED above the button flashing rapidly when you press the button for operating the heating system.

11. Reset


If a previously saved position is no longer properly set, you may have to reset the software of the Treatment bed SoSoft. To correct the differences of position, you must carry out the following steps.

NOTICE! A reset process is only possible in the event of a system error.

PLEASE PROCEED AS FOLLOWS:

Reset 1 (actuators):


Press the **lower** two buttons of the hand control (s. fig. below) **simultaneously and hold them down for approx. 10 seconds** until the signal tone ends.

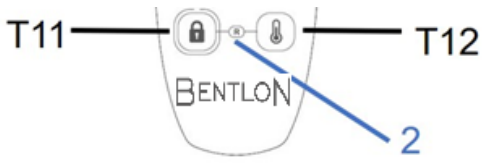

1. The backrest actuator must be completely retracted (**button T 2** of the hand control).
2. The actuator for the inclination of the seat area must be level with the leg section (calf area), i.e. set horizontally.
   ATTENTION! Please keep **button T 3** of the hand control pressed until the actuator switches off!

Reset 2 (lifting columns):

Press the **upper** two buttons of the hand control (s. fig. below) **simultaneously and hold them down for approx. 10 seconds** until the continuous tone changes to an interval tone.

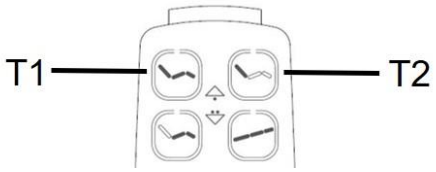

The two lifting columns must be moved to the lowest position (**button T 1** of the hand control). To do this, alternately use buttons T 1 and T 2 to lower both lifting columns steps by step until they have reached their lowest position.

- The reset process is complete when the described positions are reached.

This process must be carried out manually, i.e. each drive must be run **individually**.

.	**Please note that you keep the buttons pressed for 2-3 seconds after reaching the end position.**

⚠ **WARNING!** The adjustment (initialization) **must be** carried out manually, i.e. each drive **must be run individually**. **Do not use any of the memory buttons for this!**

⚠ **WARNING!** Ensure that there are no persons, body parts of persons, animals or objects between the lying surface and the base.
RISK OF COLLISION! RISK OF CRUSHING! RISK OF INJURY!
RISK OF DAMAGE!

12. Cleaning and Care

⚠ **WARNING** **Always** disconnect the Treatment bed SoSoft from the mains (unplug) before cleaning and maintenance. RISK OF ELECTRIC SHOCK! RISK OF FIRE!

⚠ **WARNING!** Regularly remove any dust build-up from the mains plug, the safety socket, and all electrical plug connections, e.g. with a commercially available vacuum cleaner. RISK OF FIRE!

⚠ **WARNING!** **Always** ensure that water or other liquids **cannot get** into the housing and the fittings of the Treatment bed SoSoft or into the hand control or foot controls. RISK OF ELECTRIC SHOCK! RISK OF FIRE!

NOTICE! Damage to the Spa Table, the upholstery and / or of the material of the base caused using unsuitable cleaning agents **is not covered by the guarantee or warranty!**

- **Do not wear** any arm and finger decoration such as rings or bracelets while carrying out any cleaning or maintenance operations. Otherwise you can scratch and damage the upholstery and / or the material of the base of the Treatment bed SoSoft unintentionally.
- Use water or a mild cleaning agent which is suitable for the material **only**.
- **Do not use** extremely hot water for the cleaning of the Treatment bed SoSoft.
- **Do not use** any hard sponges, scouring pads or floor clothes.
- **Do not use** any aggressive cleaning agents such as scouring powders, floor cleaners, etc.

- ❖ **Do not use** any dishwashing detergents.
- ❖ **Do not use** any solvents or bleaching agents.
- ❖ **Do not use** any essential oils.

12.1 Cleaning of the upholstery

- ❖ The upholstery of the Treatment bed SoSoft is made of **polyurethane** (**PU**).

16.1.1 Cleaning and care of the polyurethane material (PU)

- ❖ Clean the PU using a soft cloth **only** provided with a small amount of a suitable cleaning agent (also available from Bentlon: see chapter 19. UPHOLSTERY CARE AND CLEANING AGENTS).

12.1.2 Cleaning and care of the artificial leather

- ❖ Clean the artificial leather **only** with a soft cloth with a small amount of a detergent suitable for this material (also available from Bentlon: see chapter 19. UPHOLSTERY CARE AND CLEANING AGENTS).

12.2 Cleaning the housing and the fittings

- ❖ Clean the housing and the fittings using a soft damp cloth or a soft sponge and a small amount of soap only.

⚠ **WARNING! Always** ensure that water or other liquids **cannot get** into the housing and the fittings of the Treatment bed SoSoft or into the hand control or into the optional foot controls. RISK OF ELECTRIC SHOCK! RISK OF FIRE!

12.3 Cover of the upholstery

- ❖ The upholstery of the Treatment bed SoSoft is equipped with a surface made of **polyurethane** (**PU**).

- ℹ To protect these material (polyurethane) as well as the actual upholstery, we recommend the use of suitable PROTECTIVE COVERS.

PROTECTIVE COVERS AVAILABLE AT B INDUSTRIES B.V.	
Standard cosmetic applications (e.g. facials, etc.):	➢ 6 or 7 -piece terry cover.

ℹ **For detailed explanations of the protective covers see chapter 18.5 PROTECTIVE COVERS**

13. Malfunctions (possible causes and troubleshooting)

MALFUNCTION	
POSSIBLE CAUSE(S)	**TROUBLESHOOTING**
ALL MOTORS DO NOT WORK.	
1. No connection to the mains.	**1.** Check the connection status.
2. The mains supply cable is damaged.	**2.** Mains cable must be replaced (please contact the customer service).
3. Hand control / actuator plug are not plugged in properly.	**3.** Check proper connection of Hand control / plug.
4. Hand control is defective.	**4.** Hand control must be replaced (please contact the customer service).
5. Control box total loss.	**5.** Control box must be replaced (please contact the customer service).
SINGLE MOTOR DOES NOT WORK.	
1. Actuator plug is not connected properly.	**1.** Connect plug properly.
2. Hand control is defective.	**2.** Hand control must be replaced (please contact the customer service).
3. Control box is defective.	**3.** Control box must be replaced (please contact the customer service).
4. Actuator plug or plug is defective.	**4.** Cable / plug must be replaced (please contact the customer service).
5. Actuator is defective.	**5.** Actuator must be replaced (please contact the customer service).

If the Spa Table shows a malfunction not noted on the schedule above or if you have checked all items on the above schedule but the problem still exists, please contact the customer service (s. also chapter 17.1 REPAIR).

13.1 Repair

YOU CAN CONTACT THE BENTLON CUSTOMER SERVICE USING THE FOLLOWING CONTACT INFORMATION:	
CUSTOMER SERVICE:	+31 (0)495 745 005
E-MAIL:	service@bentlon.com

1. Always provide the **REF and CE number** when contacting Bentlon customer service.
2. Repairs and modifications must be carried out by B Industries B.V. or a person authorized by B Industries B.V..
3. The electrical installation of the room the Treatment bed SoSoft is installed in must meet the current legal regulations at the site where the Treatment bed SoSoft is operated.
4. The Treatment bed SoSoft is operated as described in the instruction manual.
5. In case of improper modifications or manipulations to the Treatment bed SoSoft, **any** warranty on behalf of B Industries B.V. is excluded.
6. The spa table can remain connected to the mains and does not need to be disconnected from the mains for daily use.
7. Regularly check the mains cable for signs of damage.
8. If the mains cable is damaged, it may only be replaced by B Industries B.V., or a customer service authorized by B Industries B.V. or a qualified specialist (elec-trical engineer).

⚠ DANGER! In case of defects or malfunctions of the spa table, contact customer service **immediately**. Use of the spa table **must** be dis-continued until the problem has been solved to avoid further damage to property and persons.

.	**If you wish to contact the Bentlon customer service regarding your Spa Table, we recommend taking a photo of the rating plate previously (e.g. using a mobile phone or a digital camera).** **This makes it easier for you to identify the information of the rating plate (e.g. the serial number of the Spa Table) to provide the required information to the customer service.**

14. Optionally available accessories

★ **For detailed information regarding the optionally available accessories (s. below) please contact our sales department. Our team is happy to help!**

14.1 Protective covers made of terry

- PROTECTIVE COVERS MADE OF TERRY for our Spa Tables are optionally available for an additional charge (s. fig. below).
- The PROTECTIVE COVERS MADE OF TERRY are available as a 6 or 7-piece set; with or without oil-repellent coating made of polyurethane (PU).
- The PROTECTIVE COVERS MADE OF TERRY are suitable for standard cosmetic applications.
- Our PROTECTIVE COVERS, which are specially-manufactured for our Spa Tables, will perfectly protect your working equipment from dirt and heavy use.
- Machine laundry (max. 60 °C) for hygienic cleaning is possible due to the high-quality material.
- Internal elastic bands keep the single sections exactly in the right position. All functions of the Spa Tables will continue working.
- Our COVERS are available in three different colours: white / choco and crème.
- You can order the relevant cover when purchasing the Spa Table.

★ For detailed information regarding the right size and suitable version of the covers as well as the compatibility with your Spa Table please contact our sales department. Our team is happy to help!

15. Upholstery care and cleaning agents

We recommend to clean the upholstery with water and soap only.
Cleaning products containing alcohol, solvent, bleach or abrasive products can damage this product.

16. Conditions of transport, storage, and operation

16.1 Transport

- During transport, lift the Treatment bed SoSoft using the **carrying aids only** (s. chapter 6.2.1 USE OF THE CARRYING AIDS).
- **Never** carry the Treatment bed SoSoft by its lying / sitting surface or its armrests during transport; otherwise you can damage the Treatment bed SoSoft seriously.
- Always consider the total weight of the Treatment bed SoSoft; you will need the help of several persons.
- **Always** consider the dimensions of the Treatment bed SoSoft. Narrow hallways, doors or staircases can make the transport difficult or impossible.

16.2 Ambient conditions for storage and transport

TEMPERATURE	from **-**10 °C to **+**50 °C ambient temperature

ATMOSPHERIC PRESSURE	from 700 to 1060 hPa
RELATIVE HUMIDITY	from 20% to 90% at +30 °C non-condensing

16.3 Ambient conditions for operation

TEMPERATURE	from +5 °C to +40 °C ambient temperature
ATMOSPHERIC PRESSURE	from 700 to 1060 hPa
RELATIVE HUMIDITY	from 30% to 75% at +30 °C non-condensing

17. Technical Data

GENERAL TECHNICAL DATA	
Total length	200 cm
Total width	79 cm
Weight	approx. 120 kg - 150 kg (depending on the design of the base)
Adjustable height Treatment bed SoSoft	approx. 61 cm - 91 cm
Inclination of the headrest	approx. +35°/ - 48°
Inclination of the backrest	approx. 0° - 70°
Inclination of the Trendelenburg / Anti-Trendelenburg function	approx. 14°
Lifting capacity	300 kg
Max. load capacity of the lying surface	350 kg
Max. load capacity of the armrests	25 kg

❖ Technical modifications reserved.

DRIVES INCL. CONTROLS	
Mains supply	100-240 V AC , 50-60 Hz
Power consumption	Output: approx. 350 W Standby power: approx. 0,5 W
Output voltage	24 V DC
Duty cycle (actuator)	10 % max. 2 min. continuous operation / 18 min. break
Device protection class	I
Protection class against electric shock	type B
Protection class against penetration of substances	IP X6

❖ Technical modifications reserved

HEATING SYSTEM	
Mains supply	100-240 V AC , 50-60 Hz
Power consumption	Output: approx. 120 W

❖ Technical modifications reserved.

18. Disposal

1. While unpacking, the packaging materials should be kept away from children, pets, and persons with reduced mental capabilities. RISK OF CHOKING! RISK OF INJURY!
2. At the end of its service life, this product must not be disposed with domestic waste but must be taken to a collection point for the recycling of electronic and electrical equipment. The corresponding symbol on the product must not be removed.
3. **Never** expose to open flames.
4. The Spa Table should be cleaned thoroughly before disposal to protect persons, environment, and materials.
5. Before disposing of a spa table, it is essential to **render it inoperable!** Always unplug the spa table from the mains and cut off the mains cable. Finally, you must dispose of it together with the mains plug.
6. The Treatment bed SoSoft is made of metal and plastic components, electronic components, electric cables. Dispose of the separate components in accordance with the current local disposal regulations.
7. Dispose of any paper and cardboard into the corresponding containers.
8. Dispose of any plastics into the corresponding containers
9. If suitable containers are not available at your residential area, dispose of these materials at a suitable municipal collection point for waste recycling.
10. Receive more detailed information from your retailer, municipal facilities or:

B Industries B.V.

- Phone: +31 (0)495 745 005
- E-mail: service@bentlon.com



Materials marked with this symbol are recyclable.



You can obtain all available information on waste separation from your local authorities.

19. Warranty

The following warranty conditions apply to <u>new goods</u> worldwide:

B Industries B.V. grants on all spa tables of the Treatment bed SoSoft series:

→ 10 years warranty on metal frames
→ 5 years warranty on actuators and
→ 3 years warranty on other components.

19.1 Exclusions from warranty

On delivery, it must be checked that there is no transport damage and that the accessories are complete and in proper condition. Resulting warranty claims must be registered within 8 days after delivery of the Bentlon product.

Excluded from warranty are wearing parts (such as all kinds of lamps as well as glass panels, consumables such as terry towels etc.).

<u>Furthermore; no warranty claims will be granted if:</u>

→ the spa table and its associated components have been damaged due to improper use.
→ damage to the spa table and its associated components has been caused by gross negligence.
→ the malfunction occurs through the fault of the operator.
→ stickers (such as rating plate, warning notices etc.) have been removed.
→ associated equipment or parts of the spa table or the associated equipment are opened without this being expressly permitted in the instruction manual (e.g. for cleaning purposes).
→ the instructions in this instruction manual have not been followed.
→ the spa table shows changes due to repair work by personnel not authorized by B Industries B.V..

20. EC-Declaration of Conformity

EC-Declaration of Conformity

We,

B Industries B.V.
De Noord 27
6001 DA Weert / The Netherlands

hereby declare, in sole responsibility, that, due to its design and construction, the device named below complies with the relevant, fundamental safety and health requirements of the EU.

(Type or model, name)

CE	**Lifting column beds** **Treatment bed SoSoft Bronze, Silver, Gold, Platinum**	**2023**

Relevant EU directives and harmonized standards used:

2006 / 42 / EU Machinery Directive

2014 / 35 / EU Low Voltage Directive

2014 / 30 / EU Electromagnetic Compatibility Directive

2011 / 65 / EU RoHS Directive

DIN EN ISO 12100 (2011)
Safety of machinery - general principles for design
Risk assessment and risk reduction

Responsible for documents:	Manufacturer's signature:	Place and date:
Mohammad Abbassi Gutenbergstraße 40 47443 Moers / Germany	Sammy Gharieni, CEO 	Weert, 01.04.2023

BENTLON

NOTES

THE OWNERSHIP RIGHTS TO THIS INSTRUCTION MANUAL ARE RESERVED. REPRODUCTION OR TRANSFER TO THIRD PARTIES IS PROHIBITED WITHOUT THE CONSENT OF B INDUSTRIES B.V.. B INDUSTRIES B.V. RESERVES THE RIGHT TO PROTECT ITS PROPERTY RIGHTS AND TO PROSECUTE ANY MISUSE.

SALES AND SHOWROOM	PRODUCTION, PICK-UP AND REPAIR
B Industries B.V. De Noord 27 6001 DA Weert / The Netherlands	B Industries B.V. De Noord 27 6001 DA Weert / The Netherlands